Restee PRO

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THANK YOU

Thank you for choosing our product. Please read this manual carefully before using your massage gun, and retain it for future reference. Revolutionary massage gun helps you to accelerate muscle recovery, alleviate stress and tension, and provide relief from muscle pain.

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USER GUIDE LANGUAGES

EN

To view this manual online scan this QR code or visit: www.restee.eu

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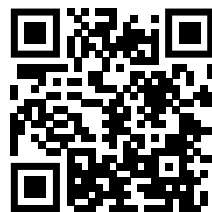
www.restee.eu

IMPORTANT INFORMATION

Incorrect operation may cause failure, electric shock, injury or fire. Before using this product, please read this manual carefully.

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PRODUCT SPECIFICATIONS

Restee PRO - Revolutionary massage gun

4x attachment heads

Battery duration: 8 Hours

Charging port: USB Type-C

Input voltage: DC 5V/1A

Power supply: lithium-ion battery

Battery capacity: 4500mAh

Recharging time: 2 hours

Item Dimensions: 13,8x8,4x4,5 cm

Weight: 450g

Material: Aircraft Grade Aluminum

SAFETY AND WARNINGS

Under certain circumstances, muscle pain may be a sign of a significant injury or illness and may require supervised medical attention and diagnosis. If you experience muscle pain lasting more than one week or if you have been in an accident, you should consult using your massage device with a medical professional.

Restee product is not intended for unsupervised usage by individuals under the age of 16. Please keep Restee product out of reach of children. Do not let children play with it.

Please note that our products are not intended for the diagnosis, treatment, or prevention of any illness or disease.

If you are pregnant or have any underlying medical conditions, you must consult with your doctor before using this product.

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CAUTIONS

Avoid immersion in water, and keep the device away from liquids and heat sources.

Use only attachment heads provided by Restee.

When the charging cable is damaged, it must be replaced with a genuine one.

Do not charge the product, in humid environments.

This appliance contains a built-in lithium-ion battery. Do not disassemble, impact, crush the battery, or throw it into fire. If there is severe swelling, do not continue to use.

Do not leave the battery in a high-temperature environment.

If you come into physical contact with any substance that may leak out of the battery, please rinse the contact area with plenty of water and seek medical attention immediately.

When charging, make sure to use a charger of 5V/1A.

As a safety feature, the device will automatically power down after 15 minutes of continuous use.

Ensure a 30-minute rest period before using the device again.

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ENVIRONMENTAL PROTECTION

Do not discard this product with general household waste at the end of its service life.

Put the product in the recycling station designated by the government.

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WARRANTY NOTICE

As a Restee consumer, you are entitled to a statutory warranty subject to the Consumer Protection Act. Further information on the consumer guarantee can be found on Restee's official website: www.restee.eu

In the unlikely event that your product has developed a fault, please contact support@resteesport.cz to arrange a replacement.

Do not attempt to disassemble any Restee devices. Warranty of your product will be voided in any instance of tampering or dismantling of the device.

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REGISTER YOUR DEVICE

To register for our Restee warranty, please fill out all of your information on our website: www.restee.eu

If anything happens to your device, please contact us with a quick photo or video of any issue you may be experiencing and your original Restee order number or proof of purchase. We will be sure to get you taken care of right away.

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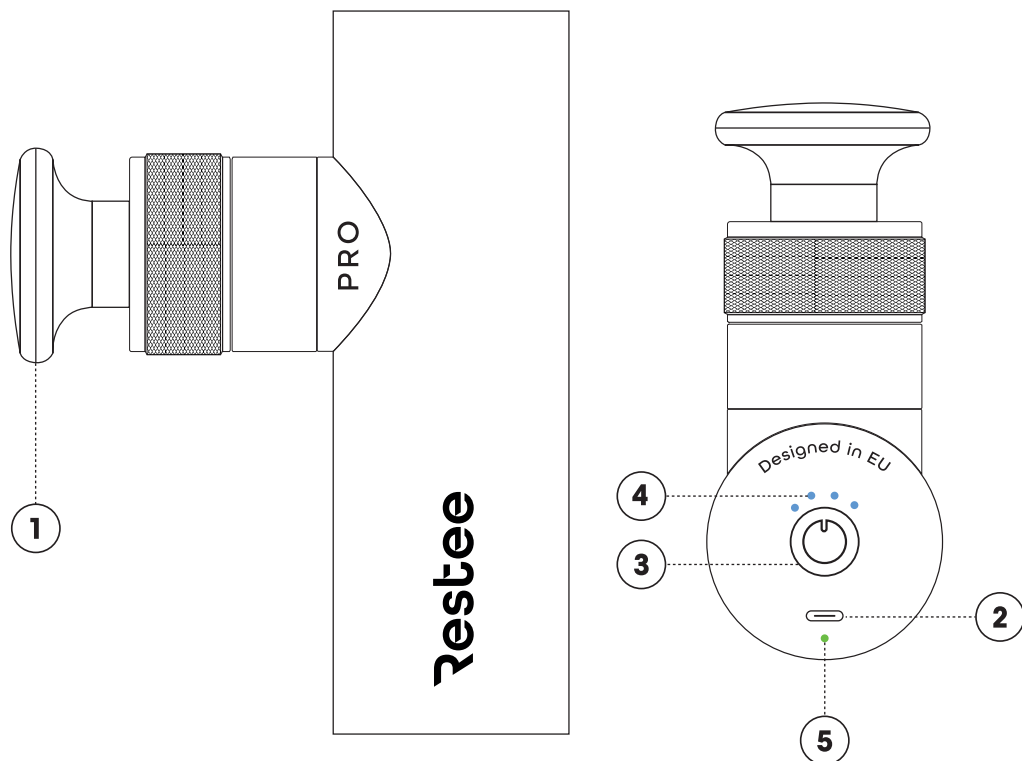
OPERATING

1. Power on device by pressing and holding the power button
2. Click the power button to select the level of speed 1-4
3. Power off the device by pressing and holding the power button

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Clean the device only with a wet cloth and dry it.

Never store the device fully charged or completely discharged for a long time. The optimal battery charge during long-term storage is around 50%.

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
CHOOSING THE BEST HEAD ATTACHMENT

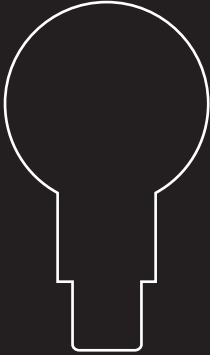ROUND HEAD

The ball-shaped head is perfectly round, allowing for smooth rotation of the point of contact with consistent and even force. It covers a larger surface area, enabling you to target the big muscle groups effectively.

Sensitive muscle groups, all areas


BULLET HEAD

Bullet head can be carried out on the soles of the feet and palm massage. Best for smaller muscle groups it allows you to dig deep into the tissues. It's great for trigger point massage but shouldn't be used for longer than 15 seconds per muscle area.

Tight muscles, deep tissue, target spots


CHOOSING THE BEST HEAD ATTACHMENT

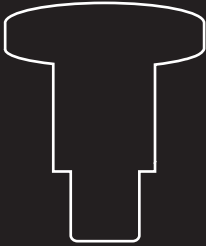FLAT HEAD

Covers a wide area for large muscle groups, while its flat top distributes the force of the massage gun. The flat round attachment is great for medium to large muscle groups. The modified force is really helpful for bulky muscle groups like the calf or biceps.

Medium to large muscles, sore muscles

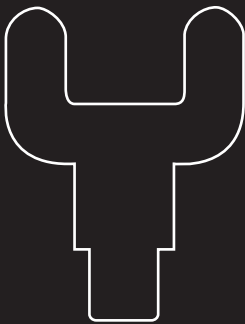
FORK HEAD

Used for trigger points and the lower back. It can be used on the sides of the spine, but always be careful to not use it directly on the bone.

Large and tight muscles, deep tissue

LEARN MORE ON OUR WEBSITE

Visit our website to discover a variety of tutorial videos. Enhance your skills and knowledge with valuable insights and tips tailored to your interests on www.restee.eu

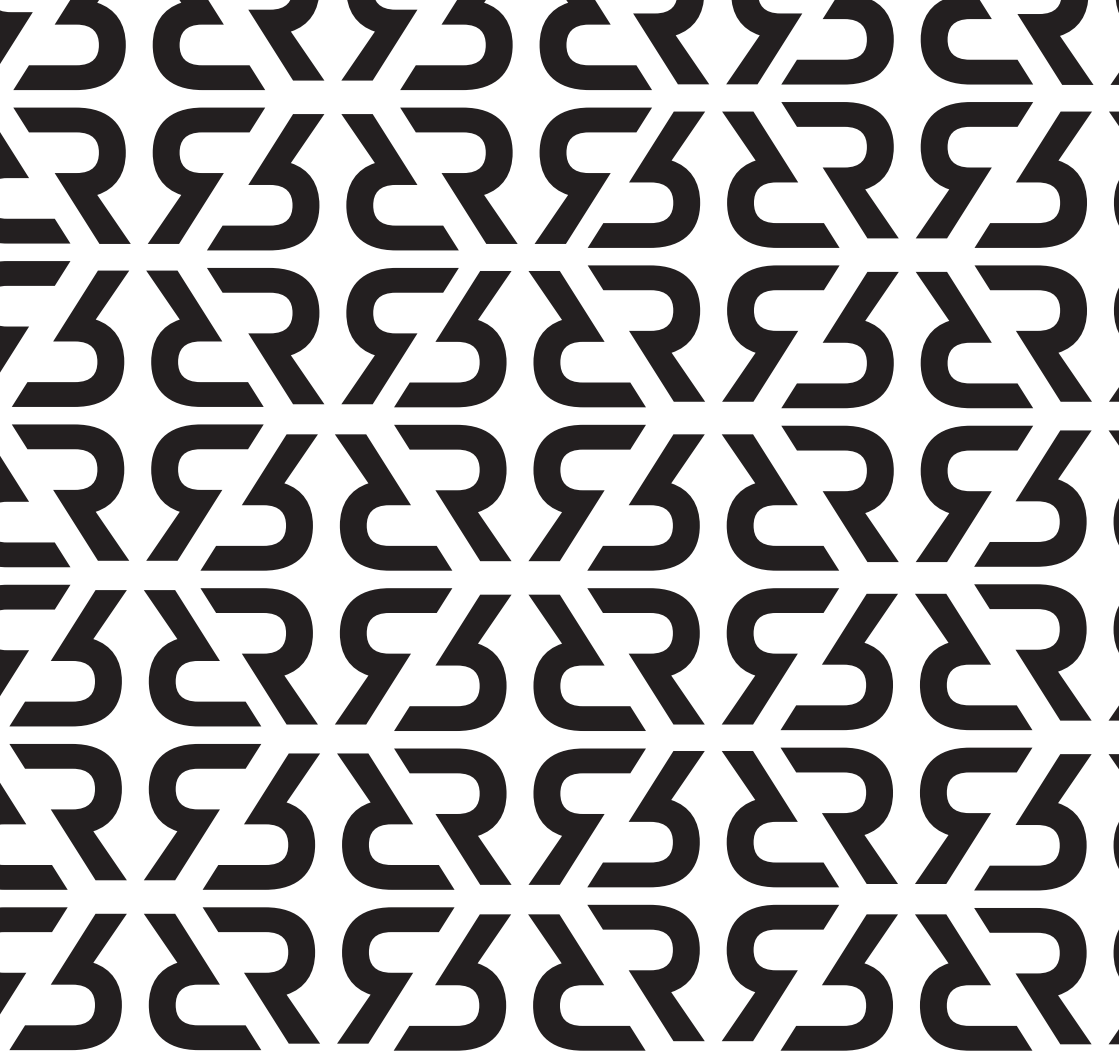Thank you for choosing Restee product.

Please read the product manual carefully before use.

For further information, please go to www.restee.eu

Socials: @restee.sport

Designed in EU by Restee

©Restee original product. All rights reserved.